

- Digital 'Pocket' Refractometer -

Improved Repeatability with

ModeS™ Technology

Pocket PAL-S

Cat.No.3860



Let's try ModeS™




Mode S™ is the solution to fluctuating measurement values.*1 Get consistent results every time!

Sample examples

Ideal for fatty, dark, and/or emulsion samples, such as dairy products.

Comparison with conventional model (PAL-1)

ModeS™ improves repeatability of samples that cause inconsistent readings.

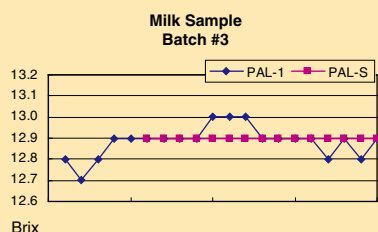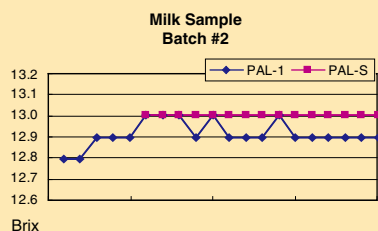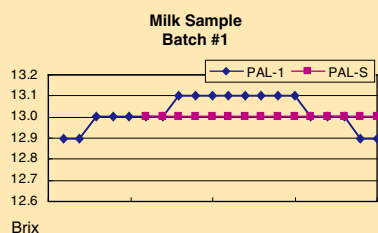

Test method : 3 batches of a milk sample were measured multiple times to test the repeatability.

Test results with PAL-1

Batch	#1	#2	#3	Overall
Average	13.0	12.9	12.9	12.9
Maximum	13.1	13.0	13.0	13.1
Minimum	12.9	12.8	12.7	12.7
Range	0.2	0.2	0.3	0.4
Standard deviation	0.07	0.08	0.08	0.10

Test results with PAL-S

Batch	#1	#2	#3	Overall
Average	13.0	13.0	12.9	13.0
Maximum	13.0	13.0	12.9	13.0
Minimum	13.0	13.0	12.9	12.9
Range	0.0	0.0	0.0	0.1
Standard deviation	0.00	0.00	0.00	0.05

Overall Results

	PAL-1	PAL-S
Average	12.9	13.0
Maximum	13.1	13.0
Minimum	12.7	12.9
Range	0.4	0.1
Standard deviation	0.10	0.05

- Conclusions -

The overall average of both units were close – 12.9 with PAL-1 and 13.0 with PAL-S.
The range of measurement values was smaller with PAL-S than PAL-1 – 0.4 with PAL-1 and 0.1 with PAL-S.
The standard deviation of PAL-S was half of that of PAL-1.



Chocolate



Ketchup



Yogurt



Caramel



Syrup



Ice cream



Mayonnaise



Blueberry jam



Condensed milk



Cutting oil



Milk



Soymilk

ATAGO®

*1: Stable readings may not be obtained with some samples due to components.

How to Measure



Place a sample on the prism.



Put the provided MAGIC™ on top.



Press the START key.

MAGIC™ (included)

This adaptor is used for a hot sample to cool it down quickly and to minimize its evaporation.



MAGIC™



The PAL takes multiple measurements for 5 to 15 seconds.
The time it takes to achieve stability varies by sample type.

Continuous
measurement for
60 seconds

Comments by our development team

It has been our mission to create a refractometer that would give good repeatability even with dark and fatty samples. In the process of developing the new algorithm to minimize fluctuations in measurement values, we conducted numerous tests with a wide variety of such samples, also taking the effects of temperature into consideration.

We are confident that the PAL-S will provide improved repeatability of samples whose measurement values may have been inconsistent with other conventional refractometers.

ModeS™ - now within your reach -

Are you tired of staring at the blurry boundary line and having to guess what the reading should be? Are you wasting time and labor on recording multiple readings and calculating the average? Try the Mode S™! The PAL-S is the first hand-held model to be equipped with this exclusive feature. The PAL-S analyzes the measurement stability internally and displays consistent readings digitally. ATAGO will continue to devote ourselves in research and development to meet the ever-changing demands of our customers.

PAL-S

Cat.No.3860

Measurement range	Brix 0.0 to 93.0% (Automatic Temperature Compensation) Temperature 9.0 to 99.9°C
Resolution	Brix 0.1% Temperature 0.1°C
Measurement accuracy	Brix $\pm 0.2\%$ Temperature $\pm 1^\circ\text{C}$
Measurement temperature	10 to 100°C
Ambient temperature	10 to 40°C
Sample	At least 0.3mℓ
Measurement time	60 seconds of continuous measurement after 5 to 15 seconds of analysis
Power supply	2 x AAA alkaline batteries
Battery life	Approximately 1,500 measurements (60 seconds of continuous measurement after 5 to 15 seconds of analysis)
International Protection class	IP65 Water Resistant
Dimensions and weight	55 (W) x 31 (D) x 109 (H) mm, 100g (main unit only)

Other models with ModeS™

Mode S™ is available on the RX-i series and selected models of the RX-α.



RX-5000 i



RX-5000 α

All ATAGO refractometers are designed and manufactured in Japan.

ATAGO CO., LTD.

Headquarters: The Front Tower Shiba Koen, 23rd Floor
2-6-3 Shiba Koen, Minato-ku, Tokyo 105-0011, Japan
TEL : 81-3-3431-1943 FAX : 81-3-3431-1945
overseas@atago.net http://www.atago.net/

ATAGO U.S.A., Inc.

TEL : 1-425-637-2107 customerservice@atago-usa.com

ATAGO INDIA Instruments Pvt., Ltd.

TEL : 91-22-2833-8038 / 8076 customerservice@atago-india.com

ATAGO (THAILAND) Co., Ltd.

TEL : 662-982-8718-9 customerservice@atago-thailand.com

ATAGO BRASIL Ltda.

TEL : 55 16 3916-6000 customerservice@atago-brasil.com

ATAGO ITALIA s.r.l.

TEL : 39 2 36557267 customerservice@atago-italia.com

ATAGO CHINA Guangzhou Co., Ltd.

TEL : 86 20-38106065 info@atago-china.com



※Specifications and appearance are subject to change without notice.

NOVATEK

ANALITİK SİSTEMLER A.Ş.

NOVATEK ANALİTİK SİSTEMLER A.Ş.
İcadiye Mh. Nevcihan Sk. No: 8/1, 34674 Üsküdar / İSTANBUL
Tel: +90 (212) 224 25 26
Faks: +90 (212) 224 00 84
info@novatekanalitik.com